

CryoStore™ Multi-Chamber Bag

The OriGen CryoStore™ Multi-Chamber Bag is indicated for freezing blood components. It provides you with the traceability, durability, and safety that you need.

Traceability

Each chamber of the bag is printed with a unique number for easy identification. The numbers on the chambers match those on the tube, assuring traceability after the aliquot process.

Durability

The spike port and film used in every CryoStore™ Multi-Chamber Bag are made from a rugged EVA. Our customers have had great success with the strength of the bag over time.

Safety

Tube sets have luer-actuated needle-free injection ports. These swabable ports help maintain sterility of the sample through numerous actuations, and help eliminate the possibility of needle sticks during processing.

Ease of Use

The standard tube set has a syringe included for removal of air prior to filling the bag with any fluid. This will help expedite the debubbling process. Each multi-chamber bag has tubing compatible for sterile welding, needle free access port and male and female luer connections.

All chambers are filled from a single location and the channel at the bottom of the bag allows to easily heat seal between each chamber with a standard tube sealer.



CS30M9




CS24M9

Technical Information


Product Code	Width (cm)	Height (cm)	Depth (cm)	Freeze/Fill Volume per Chamber	Freeze Volume	Fill volume	Max Number Chamber
CS30M9	10.8	5.5	0.9	6mL	Up to 30mL	Up to 30mL	5
CS24M9	8.4	5.5	0.9	6mL	Up to 24mL	Up to 24mL	4

CryoStore™

Multi-Chamber Bag

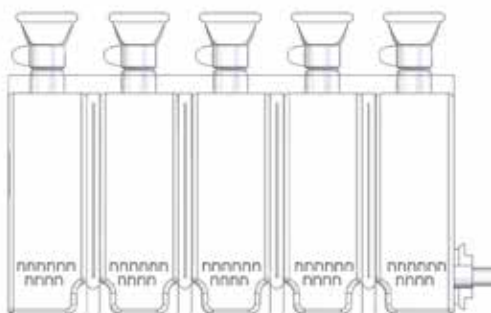

CS30

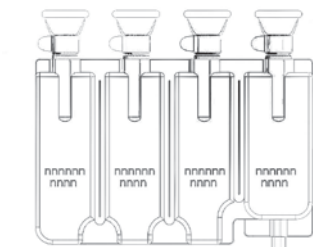

CS24



CS30M9 with empty syringe



M9 Tube Set with empty syringe

Tubing Configurations

Suffix	Tube Set Description
M9	One male luer with cap, two female luers with caps, syringe and a needle-free injection port

Packing Configurations

Product Code	Bags per Tyvek Pouch	Pouches per Case	Stock/ Non-Stock
CS30M9	1	12	Stock
CS24M9	1	12	Stock

The OriGen CryoStore™ Multi-Chamber Bag is indicated for freezing blood components. Carefully read Instructions for Use before using this product. Sterile, single use only. CryoStore™ bags are made in the USA.



7000 Burluson Rd, Bldg D
 Austin, TX USA 78744-3202
 Tel: +1 (512) 474 7278
 Tel: Toll Free (800) 233 9014
 E-mail: sales.us@origen.com