


proPAG^E




The proPAGE mini vertical system is the perfect unit for **routine vertical electrophoresis** using pre-cast or hand-cast gels. Available as a 2 or 4 gel set up for improved versatility.

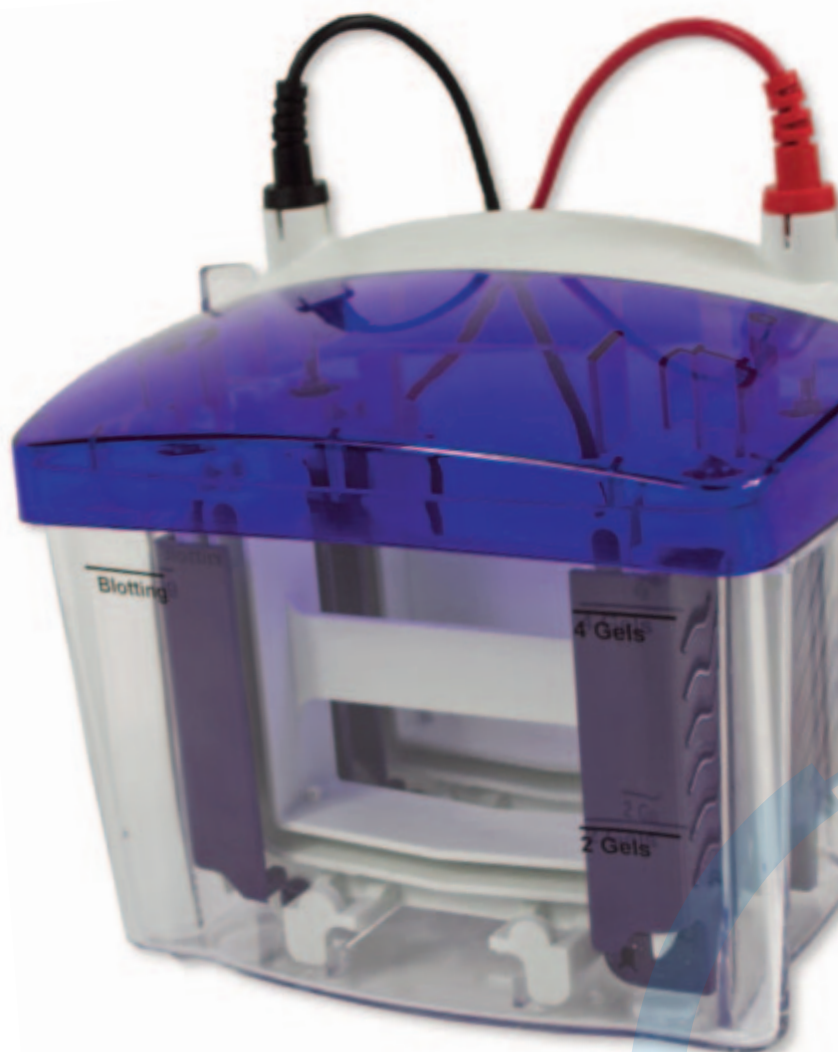

- Run 1 to 4 precast or hand-cast mini gels
- Simple and innovative Push-to-Lock cams for precision alignment of casting sandwich on any flat surface
- Bonded glass spacers available in different thickness
- Lift-up seal combined with extra soft silicon mat for leak free gel casting
- Streamline individual casters to maximise casting in restricted bench area of up to 4 gels simultaneously
- Injection Molded Mini Unit - durable, leak-proof environment for complete safety and long life

Applications examples:

- Polyacrylamide gel electrophoresis (PAGE or SDS-PAGE)
- 2-D gel electrophoresis
- Screening samples Protein profiling
- Evaluating sample preparation conditions

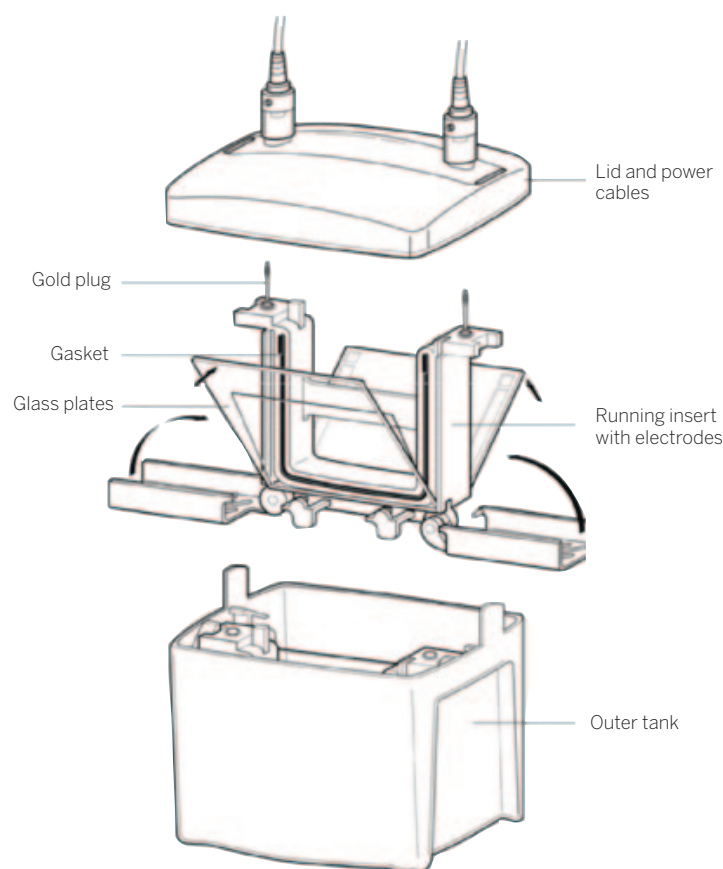

TECHNICAL SPECIFICATIONS

UNIT DIMENSION (LxWxH)	180 x 130 x 160 mm
NUMBER OF GELS	1 - 4
ACTIVE GEL DIMENSION (WxH)	83 x 73 mm
GEL THICKNESS	0.75 mm, 1.0 mm, 1.5 mm
COMBS	10 WELLS, 15 WELLS
SAMPLE CAPACITY	30 SAMPLES, 15/GEL
TANK BUFFER VOLUME	2 GELS: 700ML 4 GELS: 1000ML
PRE-CAST COMPATIBILITY	BIO-RAD, PIERCE, INVITROGEN
TYPICAL RUN TIMES FOR SDS- PAGE	35-45 MIN (AT 200 V CONSTANT)
RECOMMENDED POWER SUPPLIES	NANOPAC-300P, NANOPAC- 500 OR POWERPRO300

ORDERING INFORMATION

CVS82BSYS0.75 proPAGE 2 - mini PAGE system for 2 gels, 0.75mm
 CVS82BSYS proPAGE 2 - mini PAGE system for 2 gels, 1mm
 CVS82BSYS1.5 proPAGE 2 - mini PAGE system for 2 gels, 1.5mm

CVS84BSYS0.75 proPAGE 4 - mini PAGE system for 4 gels, 0.75mm thick
 CVS84BSYS proPAGE 4 - mini PAGE system for 4 gels, 1mm thick
 CVS84BSYS1.5 proPAGE 4 - mini PAGE system for 4 gels, 1.5mm thick

Combs

CVS8-10-0.75 Comb - 10 Sample, 0.75mm, pk/5
 CVS8-10-1.0 Comb - 10 Sample, 1mm, pk/5
 CVS8-10-1.5 Comb - 10 Sample, 1.5mm, pk/5

CVS8-15-0.75 Comb - 15 Sample, 0.75mm, pk/5
 CVS8-15-1.0 Comb - 15 Sample, 1mm, pk/5
 CVS8-15-1.5 Comb - 15 Sample, 1.5mm, pk/5

THISTLE SCIENTIFIC LTD

Unit 41, Somers Road Industrial Estate,
 Rugby, CV22 7DH
 United Kingdom

T_ +44 (0)1788 565300
 E_ ENQUIRIES@THISTLESCIENTIFIC.CO.UK
 W_ WWW.THISTLESCIENTIFIC.CO.UK